


Setting up your light

STEP ▶

1. You can choose only red or only IR or both by selecting the "Mode" button. The red indicator light will illuminate indicating which wavelength configuration is selected.
2. Adjust the time by pressing the "Time" button.
3. Press "▶||" to run the light, start the treatment.



Remember that near infrared LEDs emit light that is NOT visible to the naked eye.

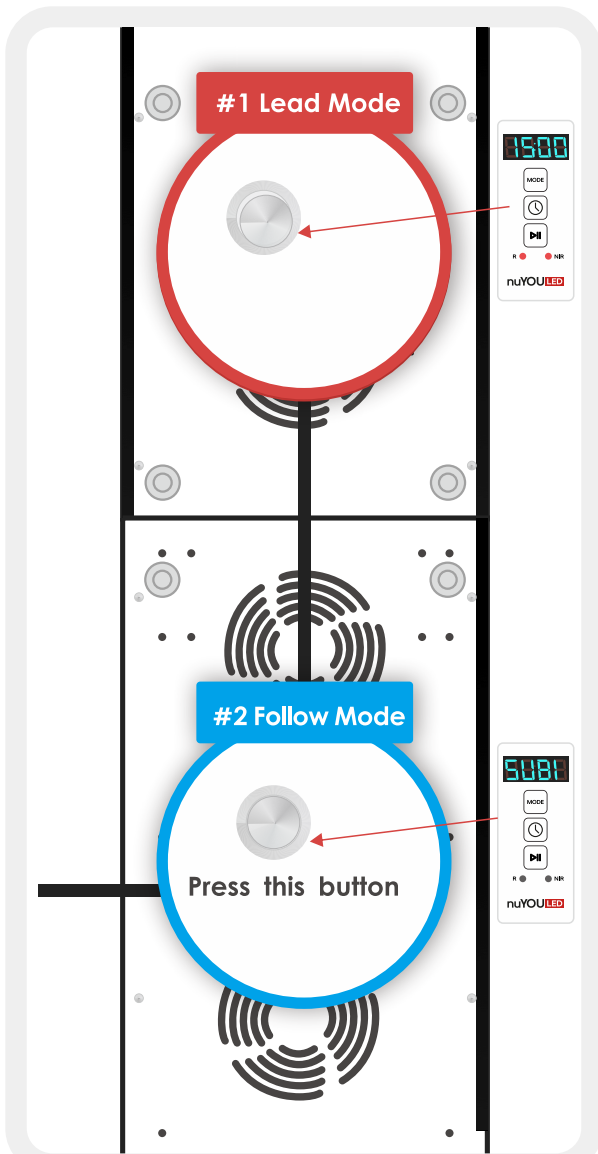

Connecting multiple lights

STEP ▶

1. Choose #1 light as leader.
2. Press #2 light's follow button, the screen will show "SUB1", #2 lights stay in "Follow Mode".
3. You only need to control the lead light, other lights will follow the lead light's instruction.