


nuYOULED

# nuMINI User Guide

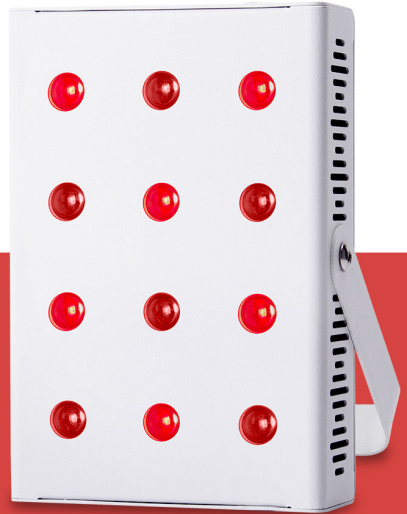

## HOW TO USE

### STEP 1

Press the power button to start treatment.

### STEP 2

The static setting when powering the nuMINI is 50/50 R/NIR.

Note: The nuMINI provides up to ten, 10 minute treatments on a full battery.

## HOW TO RECHARGE

Plug the included AC cable into the nuMINI, located on the bottom of the device and plug the other end into any standard wall outlet.

Note: The nuMINI will take 3 hours to fully charge. The nuMINI will provide 2 hours of red light treatment when fully charged.

## Treatment Instructions

2-15 minutes a day for maximum benefits. 2 to 4 days a week.

## Parts List

1 x nuMINI	1 x AC cable
1 x Travel Case	1 x Eye Protection Goggles



nuYOULED.com



@nuYOULED



@nuYOULED

LED Power Class:	60W	Cooling Fans:	n/a
Input voltage:	AC 100~240V	Frequency:	50-60Hz
Power Factor:	>0.95	Operating Temp.:	0-107F
Dimensions:	5.9" x 3.9" x 1.4"	Weight:	1.9LBS
Power Consumption:	14W	Warranty:	1 Year
LEDs:	12 x 5W	Lifespan:	50,000 hrs
Beam Angle:	30 Degree	EMF Emission:	0.0 $\mu$ T @ 6"
Best For:	Targeted Treatment	Spectrum:	Red / NIR
Certifications:	CE ROHS FCC	FDA Class II	

#### Technical & Safety Information

##### Disclaimer

nuYOU LED products are not intended to diagnose, cure, treat, or prevent any disease. Our products are general wellness/fitness and low-risk products that do not need any FDA clearance; In accordance with the "General Wellness: Policy on Low- Risk Devices" draft released on January 20, 2015.

##### Safety Considerations

Using any nuYOU therapy device against the described way can automatically disqualify you from claiming warranty. Indeed, it will result in a severe injury. nuYOU LED lights are electric devices, so for preventing electrical shock, you must follow these safety instructions:

- Do not get wet or submerge your device underwater.
- Avoid using water or any liquid while cleaning your device.
- Make sure that your hands are dry when touching the switches or plugs.
- Wipe your hands before using or touching any part of the device.
- Do not pull the cord or lift the device with the help of the cord.
- If the cord is damaged, do not turn on the device.
- Do not use the device if it has been left outside, dropped, or has any encountered with liquid or water.
- Do not use your nuYOU LED device near the heated surface.
- Do not try to repair or modify any part of the light. It does not include user-serviceable parts.
- Do not store or use your device in hot temperatures that exceed 130°F or 54°C.

#### Warranty Information & Return Policy

##### 60 Day Satisfaction Guarantee

When you choose nuYOU LED device, your satisfaction will be guaranteed. nuYOU LED gives the freedom to its buyers that they can claim their money within 60 days if they are unsatisfied for any reason.

##### 1 Year Complete Warranty

if your light is not working correctly, within one year of your purchase, we will fix or replace it. Each case will be assessed on the basis of the condition of the device. We will make sure your back up as quickly as possible.

1. If the product is non-defective, it is the responsibility of the customer to return the cost to nuYOU LED relating to shipping the product. International shipping charges are non-refundable. The customers are required to return the product in the original packaging, including all accessories. After receiving and inspecting for missing accessories or any damage to the product, nuYOU LED issues the refund. The product is sent back to the customer if the light is damaged or returned in un-sellable condition. The same payment method will be used for the refund that was used in the original order. A 20% restocking fee will be assessed on any elective returned purchases.
2. nuYOU LED will assess each warranty claim on a case by case basis. Resolution of warranty claims may include but are not limited to phone or email support to troubleshoot the product, shipment of a replacement part or parts, or replacement of the product with a new or refurbished unit. The warranty is for one year, and that will start from the date of purchase. Damage lights cannot be claimed under warranty. Any product damage or break from dropping, opening without prior authorization, and misuse will not be treated under warranty.
3. If the customer refused to accept the order that was shipped correctly, or the customer cannot receive it on his address, the order will be returned to the warehouse. nuYOU LED will charge \$50 for warehouse, shipping, and receiving charges.
4. Warranty claim is available for the original purchaser who purchased directly from [www.nuyouled.com](http://www.nuyouled.com). Authorized re-seller are non-transferable.